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Homend
Serving life

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WEBODY
User's Manual

Homend
Serving life

We know that you are eager to use your new scale very soon. However, before learning the physical values of your body, we recommend that you read the user guide we have prepared for you to get more accurate and healthy results in measurements. It will only take you 10-15 minutes to read it and then you will be able to use your scale more properly.

We have tried to make everything as clear as possible in the user guide. We hope that reading it once will be enough. We recommend that you keep this guide in a place where you can easily find it, in case something sticks in your mind in the future.

WEBODY 3903H

WHAT IS A BLUETOOTH SMART SCALE?

Bluetooth Smart Scale is a multifunctional device that not only measures your weight, but also provides different data that will allow you to optimally manage your body weight. It operates with the Bioimpedance method commonly used in home scales to calculate the weight composition of your body.

TECHNICAL FEATURES OF YOUR SMART SCALE:

Maximum carrying capacity: 180 kg

+ Automatic opening and closing functions

Note: Your smart scale will automatically turn on when a weight of 5 kg or more is placed on it.

+ Different data measurements such as body weight, body mass index, body fat percentage, muscle mass, visceral fat index, body water content, bone mass, protein content, basal metabolic rate, metabolic age, weight monitoring,

body fat content, and body weight excluding fat

+ High-precision strain gauge sensor system and high-performance processor

+ Secure and shatterproof glass push pad.

+ Warning system with low battery and overload indicators

Battery

Power source for your smart scale: AAA battery

Before first use, open the battery compartment cover of your smart scale and insert the battery, observing the correct polarity.


2. How do you use your smart scale?

Measurement of Body Weight

If you do not have a Bluetooth connection to your smart scale, only your body weight will be measured when you stand on the scale. For a full-fledged measurement, you need a Bluetooth connection. For weight measurement, follow the procedures below.

+ Place your scale on a hard, solid, and flat surface. Do not use on soft surfaces such as carpets

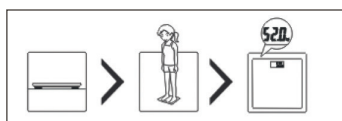+ Step onto the scale surface and remain as still as possible until your body weight appears and stabilizes on the LCD display.


+ If you take no further action, the scale will switch off automatically after a while.

Establishing the Bluetooth Connection

Download the Homend Smart app

Anyone with an iOS or Android device can download the Homend Smart app on their phone. To download the app on your phone, you must have version 8.0 for iOS phones, 5.0 for Android phones with Bluetooth 4.0 or higher.

Connect to Homend Smart App with your Smart Scale

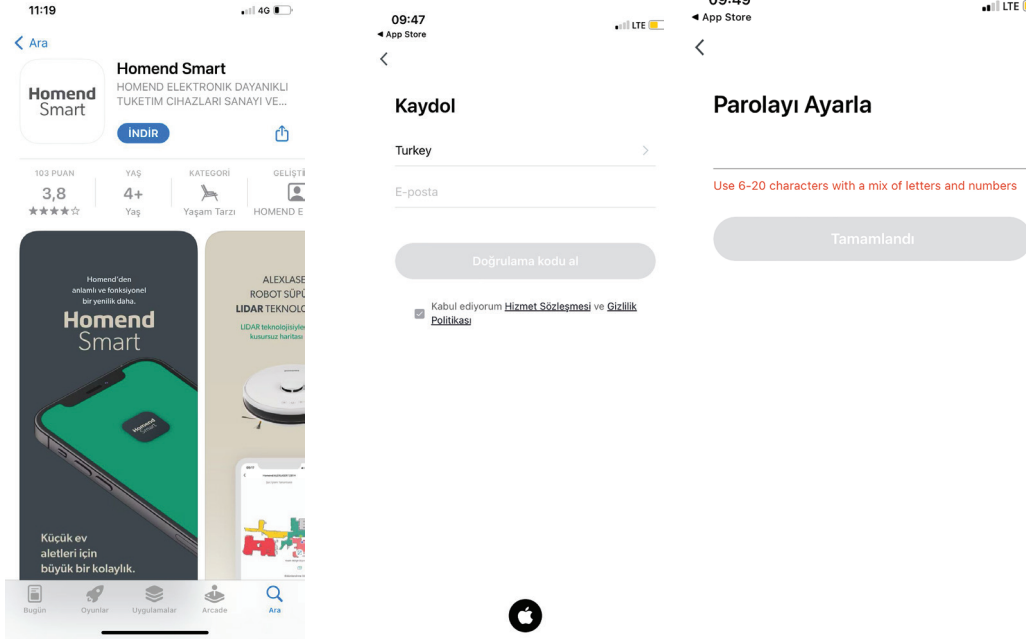
Step 1: Download the
"HOMEND SMART APP"
on your mobile device.


You can download the application
from the App Store or Google
Play.

Step 2: Register for the app and
enter the verification code that
you received via email.


Step 3: Create a password
for our app account.

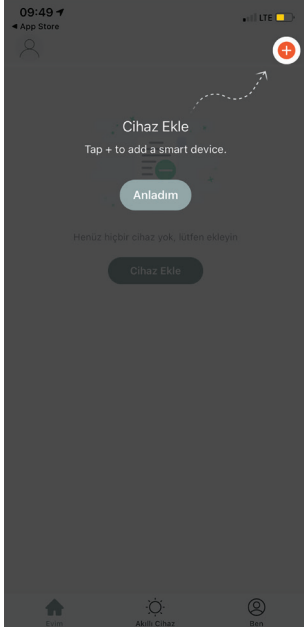
Step 4: Turn on your device.
After your device's bluetooth is activated, click on add device.

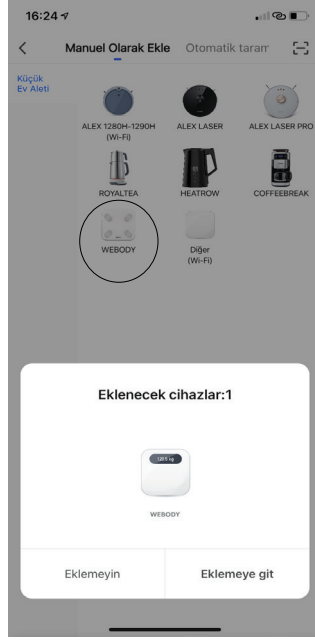
Step 5: Click on Go to add

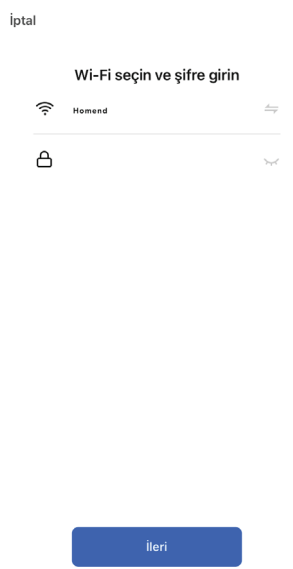Step 6: Check the name of the Wi-Fi to which your phone is connected, enter the password for your Wi-Fi network and click ileri (next).

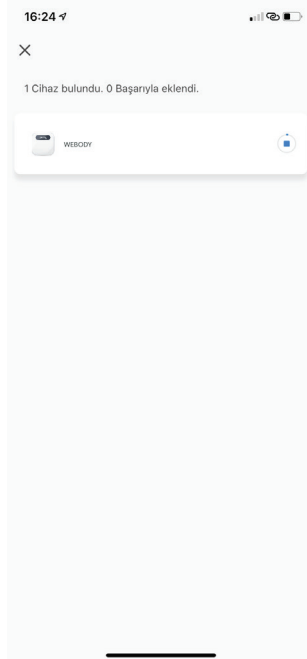Adım 7: Your device is being added, please wait

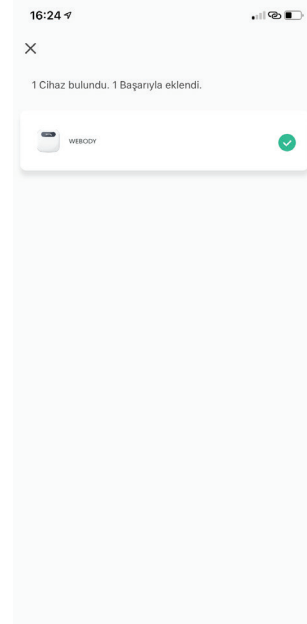Step 8: Your Smart scale has been successfully added.

Start Using Your Smart Scale with the App

When you go back to the app home page, the device number <xxx bluetooth scale> will appear on the screen. When you press this number of the device you have just paired, the (Measurements) page will open. At this time, you should stay on the scale until the measurement is complete. Once your weight measurement is fixed, you will be presented with 19 different measurement data on the screen.

The results of every measurement you take on your smart scale will be saved. You can review your measurement history on the (Records) page. The app can record and store up to 100 different data throughout the year. You can change your personal information and measurement units as you wish on the My Profile page.

You can add other users to your scale and select from existing users by clicking on the My Profile icon. You can also review your data when you are offline. However, a maximum of 20 measurement results can be stored in offline data.

When you use your smart scale without opening the app, offline data will be generated, which will then be transferred to the app.

Attention: Your scale and your phone must be paired for offline data to be transferred to the app.

Error signs that may appear on your smart scale

Overload error: When the maximum load capacity of your smart scale is exceeded, an “Err” warning will appear on the LCD display. Take care not to exceed the maximum load capacity to avoid damaging your scale.

Low battery: : If “Lo” appears on the LCD display, this indicates that the battery level in your scale is low. Open the battery compartment cover and insert a new battery into the compartment after removing the old one, observing the polarity.

Erroneous measurement:

If the error sign “Err1” has appeared on the screen of your smart scale, there are two possible reasons for this:

1. Your body fat percentage was measured below 5% or above 50%.
2. The measurement failed and could not be performed.

Do's and don'ts for proper use

- 1) For measurements of Body Weight, basal metabolic rate, muscle mass, water content, fat content, bone mass, and metabolism, you should touch the electrodes on the surface of your smart scale with bare feet.
- 2) Remember to check the batteries in case your smart scale is not running properly.
- 3) When cleaning the surface of your smart scale, use a soft cloth and alcohol or glass cleaner. Do not use soap and similar chemicals.
Do not expose your smart scale to water and do not expose it to very high or low temperatures.
- 4) Keep in mind that your smart scale is a high-precision measuring device. Never jump on the scale or disassemble the device. Handle it as carefully as possible to avoid damage when carrying it.
- 5) Your smart scale is intended for home use only and is not for commercial use. All measurements it makes are intended to give you an idea. Always consult a physician before following a diet or exercise program based on these results.

Attention: The surface of your smart scale can be slippery if it gets wet.

Model no	Webody3903H
----------	-------------

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The manufacturer's or Importer's:

Business Title: HOMEND UK LIMITED

Address: 1 Holtspur Lane, 1st Floor Wooburn Green, H.Wycombe, Buckinghamshire, United Kingdom, HP10 0AA

Customer Services: 0800 044 8181

Web: www.homend.com.tr

HOMEND UK LIMITED
Company Number 14162914
C/O Srm Business Consulting Ltd.
1 Holtspur Lane, 1st Floor Wooburn Green,
H.Wycombe, Buckinghamshire, United Kingdom
HP10 0AA

The seller's

Business Title:

Address:

Phone Number:

Authorized Person:

Seal and Signature of the Authorized Person:

Invoice Number:

Invoice Date:

Date and Place of Delivery:

The product's

Type:Weight

Model: Webody 3903H

Brand: Homend

Serial Number:

GARANTİ BELGESİ

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Unvanı: HOMEND UK LIMITED

Merkez adresi: 1 Holtspur Lane, 1st Floor Wooburn Green, H.Wycombe, Buckinghamshire,
United Kingdom, HP10 0AA

Telefonu: 0800 044 8181

Yetkili firma kaşesi / imzası:

HOMEND UK LIMITED
Company Number 14162914
C/O Srm Business Consulting Ltd.
1 Holtspur Lane, 1st Floor Wooburn Green,
H.Wycombe, Buckinghamshire, United Kingdom
HP10 0AA

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Unvanı: Karaca Porzellan Deutschland GmbH
Deutschland GmbH