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
Homend
Serving life

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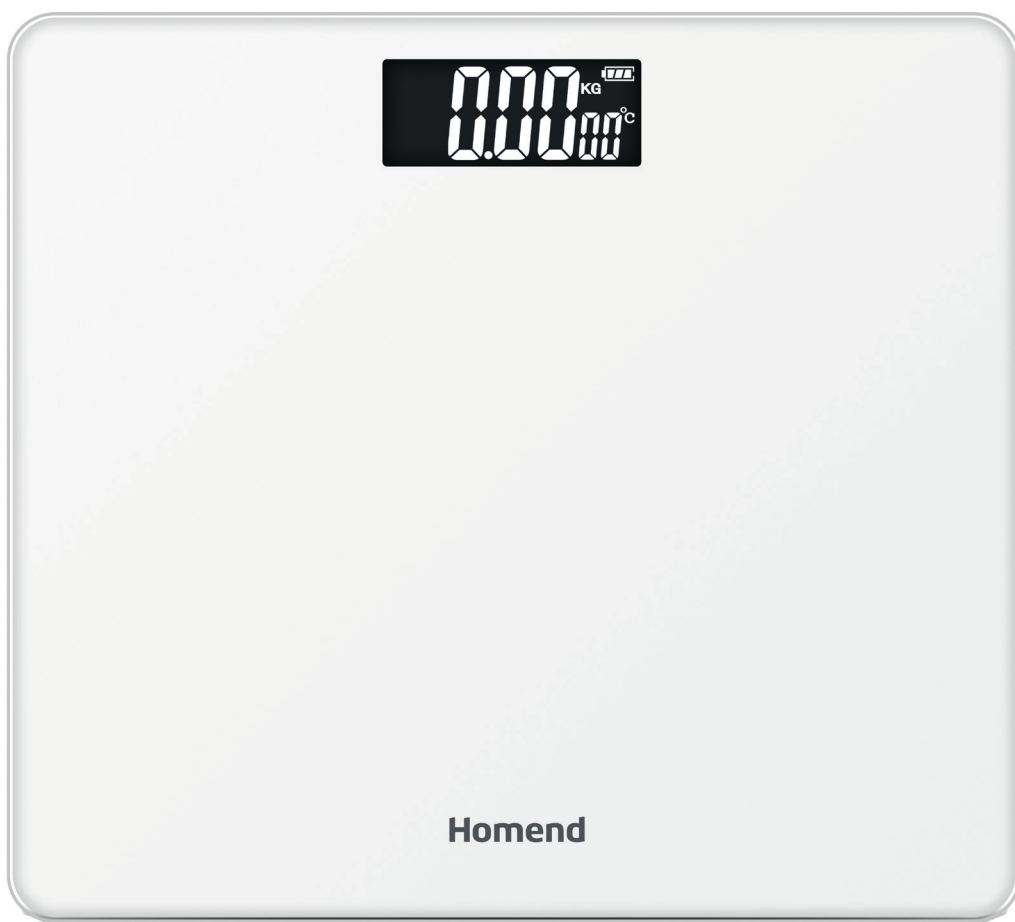
WEBODY
User's Manual

Homend
Serving life

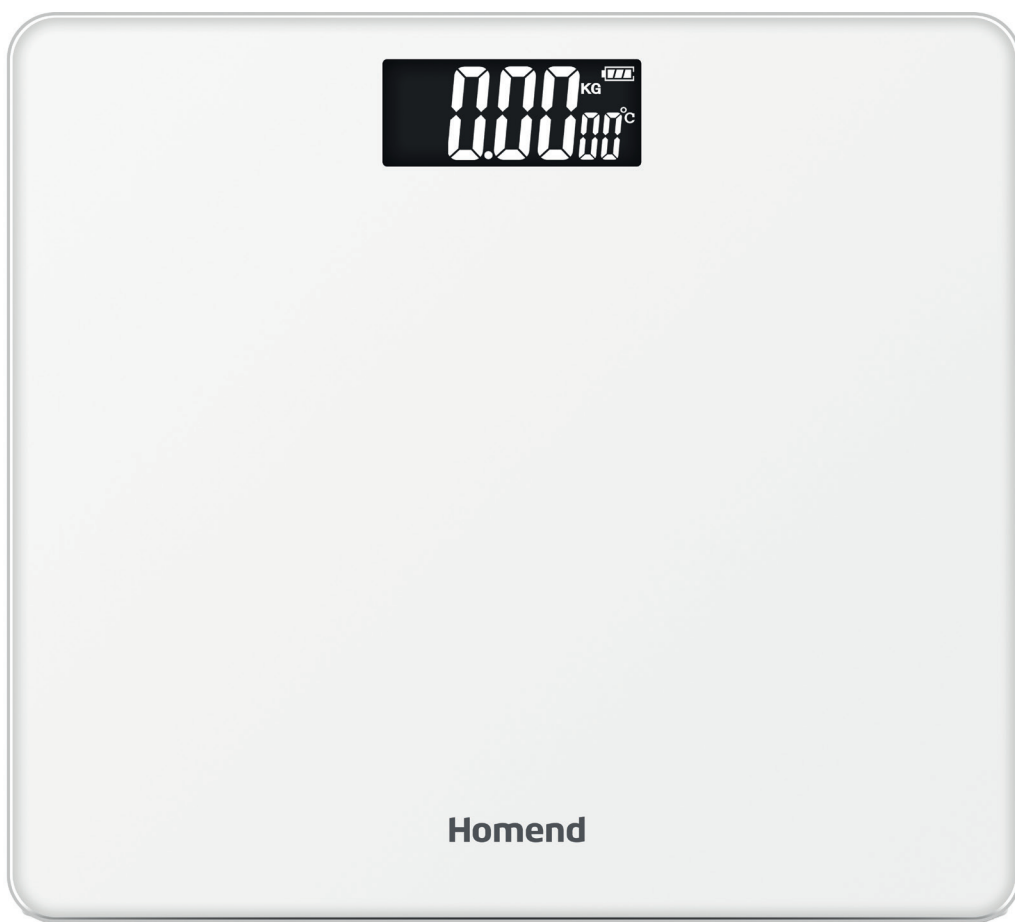We know that you are eager to use your new scale very soon. However, before learning the physical values of your body, we recommend that you read the user guide we have prepared for you to get more accurate and healthy results in measurements. It will only take you 10-15 minutes to read it and then you will be able to use your scale more properly.

We have tried to make everything as clear as possible in the user guide. We hope that reading it once will be enough. We recommend that you keep this guide in a place where you can easily find it, in case something sticks in your mind in the future.

WEBOODY

PRODUCT FEATURES

- + Display: LCD Display
- + Maximum Capacity: 180kg/396lbs
- + Sensitivity: 0.1kg/0.2lbs
- + Unit: kg/lb./st.: lb.
- + Auto Open and Auto Close
- + Low Power and Overload Indicator
- + Operating Ambient Temperature: - 100 C - 400 C
- Relative Humidity: 25% - 95% Relative Humidity
- + Storage and Transportation Temperature: - 100 C - 550 C
- Relative Humidity: 25% - 95% RH

Power supply

Battery: 2 AAA 1.5V batteries

Please install the batteries before using the scale.

Using

Step 1: Use the balance on a flat surface. Avoid using it on a carpet as this may affect accuracy.

Step 2: Step onto the platform, the body weight will be locked on the display when the weight is stabilized.

Step 3: The scale will automatically switch off after about 10 seconds of inactivity

ERROR CODES

1. Lo - When the LCD shows "Lo", please open the battery cover and replace the new battery with the correct polarity and negative.

2. Err - When the maximum capacity is exceeded, the display shows Err when the scale is overloaded (more than 180kg/396lbs).

RECOMMENDATIONS FOR USE AND MAINTENANCE

- + Please place the scale on a stable and flat surface during use to ensure the highest accuracy and lifetime.
- + Clean the balance with a wet cloth but do not allow water to seep into the unit and the body balance.
- + Please remove the battery from the scale if it will not be used for a long time.
- + Do not place anything on the scale when not in use.
- + Never use abrasive cleaners or chemical solvents.
- + Do not set the balance upright, only hold it in the appropriate flat position.
environment.

CLEANING AND MAINTENANCE

- + Please use a soft and dry cloth to clean the device and prevent water from getting inside. Never use chemical detergent.
- + Wipe the surface of the scale with a soft cloth or tissue paper. Wipe the body of the scale with a soft cloth after soaking it in water or neutral detergent, wring the cloth out, and then wipe it with a dry cloth.
- + Please do not immerse it in water or allow any liquid to enter into it when cleaning the product.
- + Please never use gasoline, volatile liquids, or strongly acidic or alkaline detergents.

MAINTENANCE AND STORAGE

- + Store the scale in a dry, dust-free, and well-ventilated place when not in use.
- + As this product is a precision instrument, please handle with care and avoid shaking, dropping or vibrating the product.
- + The product should be repaired by an authorized maintenance center.
- + To avoid damage to the device due to battery leakage, please remove the batteries when not using the product for a long period of time.
- + Battery quality depends on the battery and the product quality certificate does not include product damage caused by battery leakage.
- + Please do not place any objects on the surface when the product is not in use, to avoid damaging the sensors.
- + Please keep the product away from heat, direct sunlight, moisture, and sharp objects.


Model no	Webody3904H-05H
----------	-----------------

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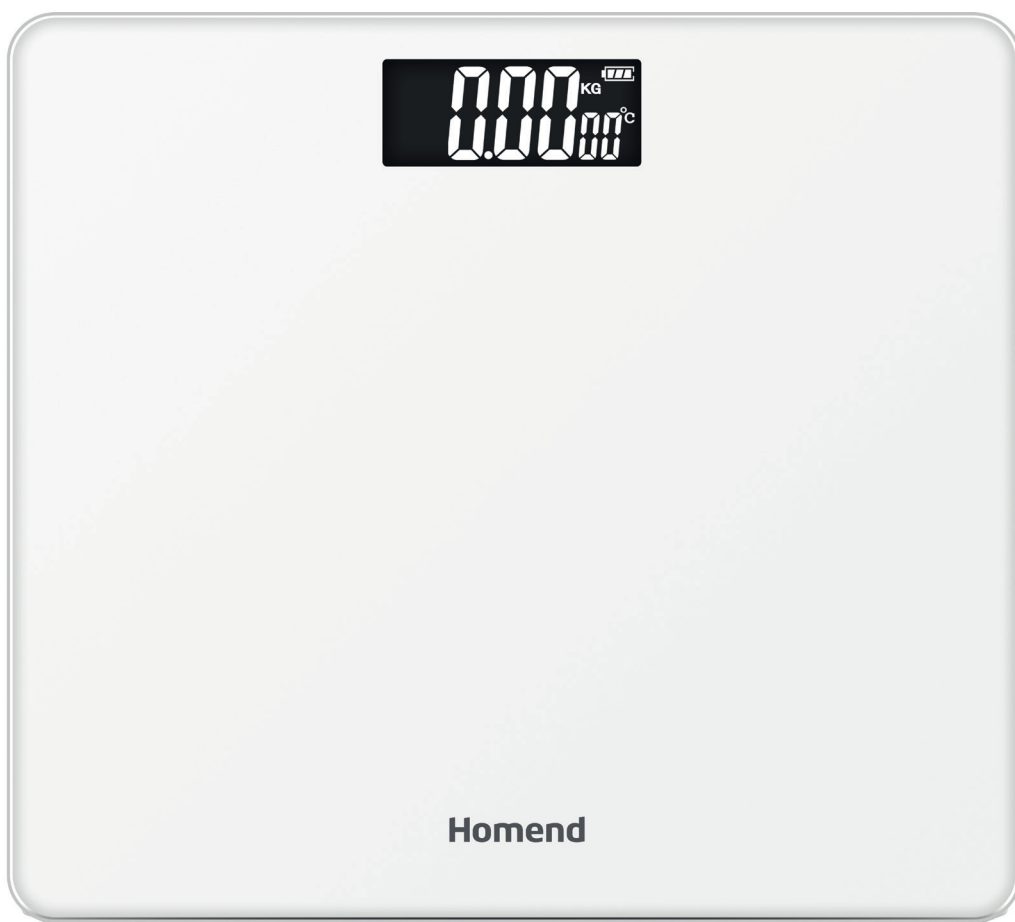The manufacturer's or Importer's:

Business Title: HOMEND UK LIMITED

Address: 1 Holtspur Lane, 1st Floor Wooburn Green, H.Wycombe, Buckinghamshire, United Kingdom, HP10 0AA

Customer Services: 0800 044 8181

Web: www.homend.com.tr

HOMEND UK LIMITED
Company Number 14162914
C/O Srm Business Consulting Ltd.
1 Holtspur Lane, 1st Floor Wooburn Green,
H.Wycombe, Buckinghamshire, United Kingdom
HP10 0AA

The seller's

Business Title:

Address:

Phone Number:

Authorized Person:

Seal and Signature of the Authorized Person:

Invoice Number:

Invoice Date:

Date and Place of Delivery:

The product's

Type:Weight

Model: Webody

Brand: Homend

Serial Number:

GARANTİ BELGESİ

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Unvanı: HOMEND UK LIMITED

Merkez adresi: 1 Holtspur Lane, 1st Floor Wooburn Green, H.Wycombe, Buckinghamshire,
United Kingdom, HP10 0AA

Telefonu: 0800 044 8181

Yetkili firma kaşesi / imzası:

HOMEND UK LIMITED
Company Number 14162914
C/O Srm Business Consulting Ltd.
1 Holtspur Lane, 1st Floor Wooburn Green,
H.Wycombe, Buckinghamshire, United Kingdom
HP10 0AA

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Kundenhotline: 0800 044 8181

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HOMEND UK LIMITED
Company Number 14162914
C/O Srm Business Consulting Ltd.
1 Holtspar Lane, 1st Floor Wooburn Green,
H.Wycombe, Buckinghamshire, United Kingdom
HP10 0AA

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Unvanı: Karaca Porzellan Deutschland GmbH
Deutschland GmbH